

The Busch Helical Method: Measuring the Ratio of an Electron's Charge to Mass

Jared Baur and Ben Sappey

(Dated: 5 March 2019)

Write abstract here.

I. OBJECTIVE

Measure the charge-to-mass ratio of the electron using the Busch Helical Method.

II. INTRODUCTION

In 1922, Busch developed a method of measuring the ratio of an electron's charge to its mass. This method involved focusing a divergent beam of electrons via a magnetic field. The electrons travel through a cathode ray tube in a specific trajectory. The trajectory of an electron in the presence of a magnetic field, when the electron is not moving parallel to the field, is in a helical motion. For mathematical purposes, this helical motion can be described as traveling along the surface of a cylinder. The electrons start from a hot cathode at the base of the tube. They pass through an anode with a velocity that depends on the potential difference between the cathode and the anode. This potential difference is called the "accelerating voltage". The applied magnetic field in the solenoid cause the beam of electrons to travel in the previously mentioned helical motion, completing a certain number of loops until hitting the fluorescent screen on the end of the cathode ray tube. The result from this process is a light that is emitted from the screen that changes with size as the magnetic field is changed. The original Busch design for this experiment had auxiliary magnets strategically placed around the cathode ray tube (instead of a solenoid) in order to produce the desired magnetic field¹.

In order to obtain the ratio of an electron's charge to its mass, let us consider a single electron moving with velocity v as it exits the anode with an angle θ to the axis of the cathode ray tube. The magnetic field will influence the electron's trajectory, however only in the motion that is perpendicular to the magnetic field. Thus, the longitudinal component of the velocity, $v \cos \theta$ will remain constant and the perpendicular component, $v \sin \theta$ will vary. The electron will go through a helical motion around the axis of the tube, all while traveling down the tube at a constant velocity. It will then pass through the center axis as it completes one helical revolution. The electron's perpendicular velocity component can be expressed as

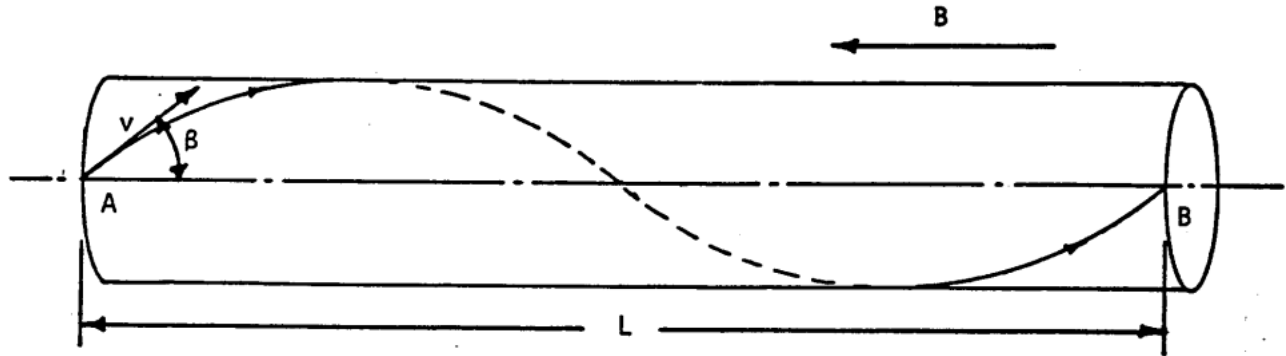


FIG. 1. Path of electrons in the solenoid. Electrons leave the cathode ray anode with a velocity v at an angle β .

$$v \sin \theta = \frac{m(v \sin \theta)^2}{eHR} \quad (1)$$

where H is the magnetic field, m is the mass of the electron, e is the charge of the electron, and R is the radius of the helical path taken by the electron. The time t for the electron to complete one revolution of this helix is expressed in Equation 2.

$$t = \frac{2\pi m}{He} \quad (2)$$

Solving for R in either of these equations and plugging into the other will give a simplified form for the time t of one revolution.

$$t = \frac{2\pi m}{He} \quad (3)$$

Since Equation 3 is independent of the perpendicular velocity, all electrons, regardless of longitudinal velocity, will complete one revolution in the same amount of time. If all electrons sent through the cathode have the same velocity (influenced by the accelerating voltage), then they will hit the fluorescent screen at the same time, thus creating a focused image on the screen. The time for the electrons to travel from the anode to the screen is expressed in Equation 4. In this equation, L is the distance from the anode to the fluorescent screen.

$$t = \frac{L}{v \cos \theta} \quad (4)$$

When the magnetic field is adjusted so that there is a focused point on the screen, the time for the electron to reach the screen and the time for it to complete one revolution are equal. Thus, Equation 3 and 4 can be set equal to one another. Solving for the ratio e/m gives equation 5. The velocity of the electron can be found by analyzing the kinetic energy $Ve = 1/2mv^2$, where V is the accelerating voltage.

$$\frac{e}{m} = \frac{8\pi V \cos^2 \theta}{H^2 L^2} \quad (5)$$

III. APPARATUS AND METHODS

to-do: name the different apparatus models

The apparatus consists of a cathode ray tube and a solenoid. The cathode ray anode emits electrons with the assistance of a high voltage power supply. The cathode ray tube is placed in the center of the solenoid with the axes of the tube and solenoid aligned. An accelerating voltage is applied to the cathode ray tube that accelerates electrons in a beam. A voltage applied to the intensity control grid reduces the current of the emitted electron beam. The solenoid provides a curved, helical path for the electron beam to follow in the cathode tube. The solenoid's current is reversible; when reversed, the magnetic field is reversed and the helical orientation of the electron beam is flipped (clockwise or counterclockwise).

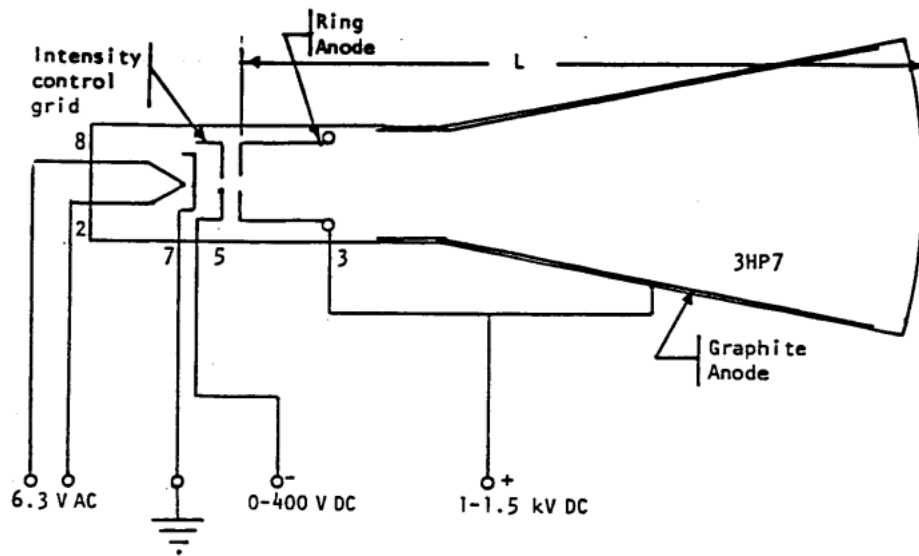


FIG. 2. Cathode ray tube schematic. Electrons are emitted from the anode and travel a distance L through the cathode tube.

Before beginning the experiment, measurements of the apparatus were made in order to obtain the values needed for calculations. The cathode ray's length L was measured (used in Equation 5) with the use of a cathetometer. The cathetometer was chosen to make this measurement since it can easily be used to measure distances of objects where the measuring tool cannot be placed close to the point of interest. The cathetometer was used to measure the distance from the anode to the end of the tube (length L shown in Figure 2). This length is defined as the distance from the anode to the phosphor screen that is the point at which the electrons are illuminated in the cathode tube. The phosphor screen is actually embedded underneath a glass of measurable thickness, which also required the use of a cathetometer to measure. In order to measure the glass thickness, the pinpoint of a thumbtack was placed on the surface of the glass. The microscope on the cathetometer was focused on this pinpoint and the reading on the cathetometer scale was taken. Next, the focus of the microscope was adjusted to the phosphor screen under the glass. Once the phosphor was in focus, the cathetometer scale reading was noted. The difference between these two readings represent the thickness of the glass, which was then used to calculate the actual length L of the cathode ray tube (from the anode to the phosphor screen).

The magnetic field in a solenoid is given by Equation 6, where H is the magnitude of the magnetic field, n is the density of turns for the length of the solenoid (number of turns per unit length), I is the current, and θ is the angle subtended at the center of the solenoid by the outer edge of the solenoid and its axis². The turn density was measured by counting the total number of turns for the entire solenoid, then dividing this count by the total length of the solenoid.

$$H = nI \cos \theta \quad (6)$$

IV. DATA ANALYSIS

V. CONCLUSION

¹Stranathan, *The Particles of Modern Physics* (Blakiston, 1942).

²*The Measurement of e/m by the Busch Helical Method*, Occidental College Physics Department (2018).

³Tipler and Llewellyn, *Modern Physics* (Worth Publishers, 1978).

⁴Taylor, *An Introduction to Error Analysis: The Study of Uncertainties in Physical Measurements*, 2nd Ed. (University Science Books, 1996).

⁵Hoag and Korff, *Electron and Nuclear Physics* (Van Nostrand, 1948).